



Hospital bed

S 962-2

Völker. Caring for People

VÖLKER

Homelike comfort

For increasing care requirements.

The S 962-2 Völker hospital bed meets a wide variety of requirements in hospitals. It enables patients to go through the various stages of mobilisation in a safe way, with support from the care staff – from sitting in the comfort chair position to getting up and leaving the bed. Also available in version S 962-2 W, which is suitable for washing tunnels.

Through a successful combination of form and function, Völker hospital beds provide greater patient comfort and enable care staff to work safely and efficiently. Völker's unique healthcare design ensures that beds fulfil their medical purpose and at the same time create a homelike ambience that supports patients' well-being and recovery.

Three convincing elements

- ✓ Split, height-adjustable side rails – move along with the mattress as it is adjusted.
- ✓ With innovative Völker MiS® Micro-Stimulation System. For improved patient comfort and better sleep.
- ✓ Also available as a wash tunnel suitable model.



**SAFE
WORKINGLOAD**



**WIDE
ADJUSTMENT RANGE**



**DYNAMIC
EXTENSION**



**LOW
BED HEIGHT**








Safe mobilisation for patients


Relief for the care staff thanks to the patented Völker side rail.

Independent mobilisation for patients helps with recovery following a period of illness. This mobilisation is supported by the central gap in the Völker side rail, alleviating the restricted freedom of movement imposed by side guards. The stable, multi-telescopic design of the Völker side rail enables patients to stand up independently and regain confidence in their own bodies. As standard, the Völker side rails are fitted on all sides of the lying surface and offer optimum protection in every position. A fixed split version of the side rail is also available as an option.

Technical data

Safe working load		250 kg
Unladen weight		approx. 143 kg
Lying surface width		90 / 100 cm
Lying surface length		200 / 210 cm
External dimensions		99 – 104 cm*

*for fixed split side rails

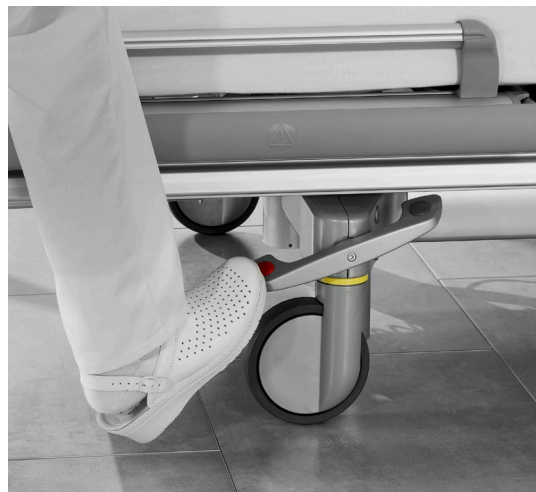

Central brake

Can also be operated at the lowest bed position.
This means that the bed can be moved at any height setting.



Learn more about
the technical data
voelker.de



Greater comfort and relaxation

The micro-stimulation lying surface.

All of Völker's lying surfaces feature the patented MiS® Micro-Stimulation System, which is based on neuroscientific principles. It reacts to even the most minute of the patient's own movements, relaying them back to the body as tactile stimuli from the lying surface through the mattress. Experience in recent years has shown that such feedback can lead to improved perception of the body, as well as improvements in sleeping and lying comfort. The combination of lying surface and mattress may also reduce the risk of bedsores.¹



Optional transverse hand control

Only available with moving split side rail.



Automatic features

Defined positions can be set quickly and easily with a simple double click. Enables users to better focus on care activities.



People requiring care have a particular right to respect and regard for their dignity. Caregivers also have a right to equipment that makes their daily work as easy as possible. We think, develop and manufacture with this in mind.

24/7

**ALWAYS
AVAILABLE**



**FOCUS
ON SERVICE**



**HIGH-QUALITY
MATERIALS**



**CERTIFIED
SUSTAINABILITY**

Head office

Völker GmbH

Wullener Feld 79
58454 Witten
Germany
+49 (0) 2302 960 960

Völker Schweiz GmbH

Moosstrasse 2
6003 Luzern
Switzerland
+41 (0) 44 563 1939

Völker Medizinprodukte GmbH

Kaufmannstraße 14
6020 Innsbruck
Austria
+43 (0) 1 205 8365

Voelker France SAS

27 Rue Maurice Flandin
69444 Lyon Cedex 03
France
+33 (0) 18 288 1530

Voelker Nederland B.V.

Laan 21
1741 EA Schagen
Netherlands
+31 (0) 20 808 1603

General enquiries

info@voelker.de
+49 (0) 2302 960 960

Sales / Service / Replacement parts

kundenservice@voelker.de
+49 (0) 2302 960 9630